

## PROVIDE RESULTS. DRIVE PROFITABILITY.

Tighten skin\* | Improve blemishes | Revitalize skin

### FACE

Before & after 3 treatments  
Courtesy of J.H. Yoon, MD, Korea,  
YonseiJams clinic



### NECK

Before & after 1 treatment  
Courtesy of S. Doherty, MD



### ABDOMEN

Before & after 2 treatments  
Courtesy of S. Doherty, MD



For more information, visit  
[www.cynosure.com/potenza](http://www.cynosure.com/potenza)

Potenza Specs & Accessories	
RF Type	Bipolar, Monopolar
Power	Max 50 watts
Frequency	1MHz, 2MHz
Repetition Rate	0.5-3.0 sec (0.1 sec step)
Handpieces	100-7043-001: Motorized Handpiece 100-7043-002: AC Handpiece
Tips	<b>Motorized Handpiece</b> 100-7043-051: S-16 100-7043-052: S-25 100-7043-053: S-49 100-7043-054: I-16 100-7043-055: I-25 100-7043-056: I-49 <b>AC Handpiece</b> 100-7043-060: P1-08 100-7043-061: A1-12 100-7043-062: A1-15
Needle Thickness	0.25
Treatment Duration	5-500ms (5ms step)
Electrical Power	100-240V, 50/60Hz, 500VA
Dimensions	19.82 x 14.40 x 14.22 inches (Incl. cable stand: 21.66 x 27.60 x 14.41in)
Weight	Approx. 21lbs
Cart	18.82 x 18.98 x 38.08 in, 71.13lbs

\*Through soft tissue coagulation.

©2020 Cynosure, LLC. All rights reserved. Cynosure is a registered trademark of Cynosure, LLC. Tiger Tip is a trademark of Cynosure, LLC. Potenza is a trademark of Jettys Medical, Inc. Cynosure, LLC. owns exclusive rights to photography. Use of photography without written permission of Cynosure is prohibited.

The Potenza radiofrequency microneedling device is intended for electrocoagulation and hemostasis of soft tissue for dermatologic conditions. Potential side effects include temporary redness, temporary tingling, and burning sensation while receiving treatment. Like all medical procedures, not all patients are suitable for the treatment. A qualified practitioner is solely responsible for evaluating each subject's suitability to undergo treatment and for informing those being treated about any risks involved with the treatment, pre- and postoperative care, and any other relevant information. Individual results may vary and are not guaranteed.

921-7027-214 Rev. 5.07/18



Potenza™  
RF Microneedling

## The new standard in RF MICRONEEDLING

The Potenza system takes microneedling to new levels. The versatility of this system allows you to customize treatments based on your patients' specific needs – **all skin types, anywhere on the body, any time of year.**

- ✓ **World's first 4-mode RF**  
Monopolar, bipolar, 1 & 2 MHz capabilities. Provide shallow and deep treatments in one system.
- ✓ **Unique Tiger Tip™ technology**  
Treat more tissue in a single treatment.
- ✓ **Single-needle handpiece**  
Target and improve blemishes.



# THE WORLD'S FIRST 4-MODE RF MICRONEEDLING DEVICE

TREAT A LARGER VARIETY OF  
YOUR PATIENTS' CONDITIONS  
**BOTH SUPERFICIAL AND DEEP**

## THE BEST OF BOTH WORLDS

### MONOPOLAR RF

Deliver energy across a larger area of tissue for deep heating, skin tightening.\*

### BIPOLAR RF

Deliver energy in concentrated, more superficial tissue. This is ideal for skin revitalization.



*Potenza is a highly advanced device for RF microneedling. It allows me to individualize and optimize treatments for my patients who seek out our practice as having the most effective and state-of-the-art technologies.*

-Dr. Robert Weiss


# A SOLUTION FOR EVERY PATIENT

Choose the right tip for the right patient



### UNIQUE TIGER TIP™ TECHNOLOGY SEMI-INSULATED NEEDLES

- 16, 25 and 49 needle densities
- Adjustable depth 0.5mm- 4.0mm

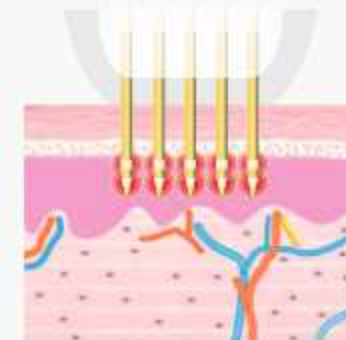


Used for treating more tissue in a single treatment without sacrificing the epidermis



### SINGLE INSULATED NEEDLES

- .08mm, .12mm and .15mm depth

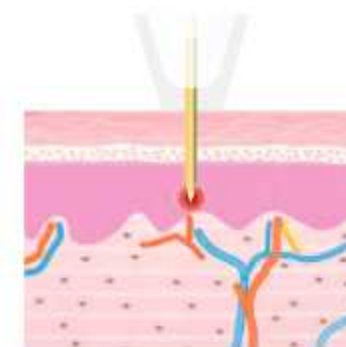


Used to precisely target and improve blemishes found in the skin



### INSULATED NEEDLES

- 16, 25 and 49 needle densities
- Adjustable depth 0.5mm- 4.0mm

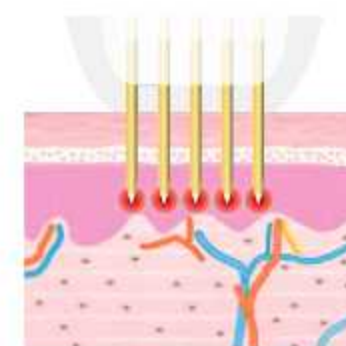

Used for both targeted and deep treatments in the skin